

# Endoscopic Cubital Tunnel

Increased safety,  
reduced recovery



## Endoscopic Cubital Tunnel (ECuTR)

SegWAY ECuTR is the first Endoscopic Cubital Tunnel Release system that allows the surgeon to safely view the entire ligament from above and below through controlled portals in the guide. It is also the only system that the scope and the blade work independently of each other giving the surgeon a better field of view and a tactile feel. The system is designed so that the Ulnar nerve may be safely viewed and followed before the release is performed. The guide covers the blade on both sides of the ligament so as to protect any unwanted release.

### FEATURES AND BENEFITS

- > **SAFETY FIRST** - **DUAL SCOPE CHANNELS** allow the surgeon to view the ligament above and below to ensure that there is no nerve
  - **SPLIT GUIDE** allows the surgeon to rotate the scope to view and follow the nerve
  - **COVERED BLADE CHANNEL** protects surrounding tissue and nerve while releasing the ligament
- > **VERSATILITY** - **ROTATING GUIDE** allows the surgeon to release with either hand or to move the release over
- > **ASSURANCE** - **LARGE FIELD OF VIEW** from above and below to make sure no nerve is in the release
- > **CONFIDENCE** - **LONGER GUIDE** allows the surgeon to release as much as necessary
- > **UNIportal** - **SMALLER** incision allows for faster recovery

Where Tradition Meets Innovative Technology




Uniportal dual scope guide system



Individual long guide system



Independent antegrade knife operation with view port

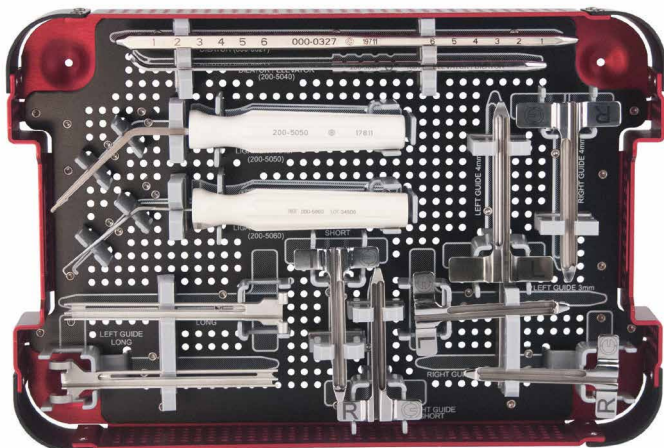

## PRODUCT LINE INFORMATION

Description	Part #
SegWAY Instrument Set (Complete Set)	200-5560
Left Medium SS Knife Guide (4mm)	200-5012
Right Medium SS Knife Guide (4mm)	200-5022
Left Small SS Knife Guide (2.7mm)	200-5014
Right Small SS Knife Guide (2.7mm)	200-5024
Left Large SS Knife Guide (4mm)	200-5010
Upper ECuTR Guide (4.0mm long)	200-6041
Lower ECuTR Guide (4.0mm long)	200-6042
Synovial Dilator / Elevator	200-5040
Rasp, Ligament	200-5050
Probe, Ligament	200-5060
Dilator	000-0327
Retrograde Ligament Knife (SINGLE USE)	200-1003
Antegrade Ligament Knife (SINGLE USE)	200-1006